

Docket No. 246424US8



IFW

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

IN RE APPLICATION OF: Thomas NILSSON, et al.

GAU: 1614

SERIAL NO: 10/728,986

EXAMINER:

FILED: December 8, 2003

FOR: A MEDICAL PRODUCT

SUBMISSION NOTICE REGARDING PRIORITY DOCUMENT(S)

COMMISSIONER FOR PATENTS
ALEXANDRIA, VIRGINIA 22313

SIR:

Certified copies of the Convention Application(s) corresponding to the above-captioned matter:

☒ are submitted herewith

☐ were filed in prior application filed

☐ were submitted to the International Bureau in PCT Application Number _____
Receipt of the certified copies by the International Bureau in a timely manner under PCT Rule 17.1(a) has been acknowledged as evidenced by the attached PCT/IB/304.

Respectfully Submitted,

OBLON, SPIVAK, McCLELLAND,
MAIER & NEUSTADT, P.C.



Richard L. Treanor

Registration No. 36,379

Fred L. Braun

Registration No. 56,123

Customer Number

22850

Tel. (703) 413-3000
Fax. (703) 413-2220
(OSMMN 11/04)

PRV

PATENT- OCH REGISTRERINGSVERKET



Intyg Certificate

Härmed intygas att bifogade kopior överensstämmer med de handlingar som ursprungligen ingivits till Patent- och registreringsverket i nedannämnda ansökan.

This is to certify that the annexed is a true copy of the documents as originally filed with the Patent- and Registration Office in a connection with the following patent application.

- (71) Sökande: Microdrug AG, 6052 Hergiswil NW, CH
Applicant (s)
- (21) Patentansökningsnummer 0303269-5
Patent application number
- (86) Ingivningsdatum 2003-12-03
Date of filing

Stockholm, 2006-05-23

För Patent- och registreringsverket
For the Patent- and Registration Office

Hjärdís Segerlund

Avgift

Fee 170:-

Ink. t. Patent- och reg.verket

2003-12-03

1

A medical product

Huvudfaxen Kassen

TECHNICAL FIELD

The present invention relates to a medical product comprising inhalable doses of tiotropium loaded in a moisture-tight, dry container and in particular, a metered dry powder medicinal dose of tiotropium bromide being adapted for administration by a dry powder inhaler device.

BACKGROUND

Asthma and chronic obstructive pulmonary disease (COPD) affect more than 30 million people in the United States. More than 100,000 deaths each year are attributable to these conditions. Obstruction to airflow through the lungs is the characteristic feature in each of these airway diseases, and the medications utilized in treatment are often similar.

Chronic obstructive pulmonary disease (COPD) is a widespread chronic lung disorder encompassing chronic bronchitis and emphysema. The causes of COPD are not fully understood. Experience shows that the most important cause of chronic bronchitis and emphysema is cigarette smoking. Air pollution and occupational exposures may also play a role, especially when combined with cigarette smoking. Heredity also causes some emphysema cases, due to alpha1 anti-trypsin deficiency.

Administration of asthma drugs by an oral inhalation route is very much in focus today, because of advantages offered like rapid and predictable onset of action, cost effectiveness and high level of comfort for the user. Dry powder inhalers (DPI) are especially interesting as an administration tool, compared to other inhalers, because of the flexibility they offer in terms of nominal dose range, i.e. the amount of active substance that can be administered in a single inhalation.

Anticholinergic agents, e.g. tiotropium, especially tiotropium bromide, are effective bronchodilators. These medicaments have a relatively fast onset and

+46 18 153050

Ink. t. Patent- och reg.verket

2003-12-03

2

Huvudfaxen Kassar

long duration of action, especially tiotropium bromide, which may be active for up to 24 hours. Anticholinergic agents reduce vagal cholinergic tone of the smooth muscle, which is the main reversible component of COPD. Anticholinergic agents have been shown to cause quite insignificant side effects in clinical testing, dryness of mouth and constipation are perhaps the most common symptoms. Because it is often very difficult to diagnose asthma and COPD correctly and since both disorders may co-exist, it is advantageous to treat patients suffering temporary or continuous bronchial obstruction resulting in dyspnoea with a small but efficient dose of a long-acting anticholinergic agent, preferably tiotropium bromide, because of the small adverse side effects.

Tiotropium bromide is the preferred anticholinergic agent because of its high potency and long duration. However, tiotropium is difficult to formulate in dry powder form to provide acceptable performance in terms of dose efficacy using prior art DPIs. Dose efficacy depends to a great deal on delivering a stable and high fine particle dose (FPD) out of the dry powder inhaler. The FPD is the respirable dose mass out of the dry powder inhaler with an aerodynamic particle size below 5 μm . Thus, when inhaling a dose of dry medication powder it is important to obtain by mass a high fine particle fraction (FPF) of particles with an aerodynamic size preferably less than 5 μm in the inspiration air. The majority of larger particles (>5 μm) does not follow the stream of air into the many bifurcations of the airways, but get stuck in the throat and upper airways, where the medicament is not giving its intended effect, but may instead be harmful to the user. It is also important to keep the dosage to the user as exact as possible and to maintain a stable efficacy over time, and that the medicament dose does not deteriorate during normal storage. For instance, Boehringer Ingelheim KG (BI) markets tiotropium bromide under the proprietary name of Spiriva®. Surprisingly, in a recent investigation into the inhalability of Spiriva® we have found that the Spiriva®/HandiHaler® system from BI for administration by inhalation of

doses contained in gelatin capsules shows poor performance and has short in-use stability.

Thus, there is a need for improvement regarding a medical product comprising inhalable dry powder doses of tiotropium bromide, for instance Spiriva®, and suitably adapted inhaler devices for the purpose of administration.

SUMMARY

The present invention discloses a medical product for use in the treatment of respiratory disorders, and comprises a metered dose of a tiotropium dry powder formulation, directly loaded and sealed into a moisture-tight, dry container acting as a dry, high barrier seal against moisture. The container itself does not emit water, which may affect the tiotropium powder inside.

Thus, the container does not release any water to the dose and ingress of moisture from the exterior into the container is thereby prevented.

The dose of tiotropium is further intended for inhalation and the container is so dry and tight that the efficacy of the dose when delivered is unaffected by moisture.

In another aspect of the invention a type of inhaler is disclosed, which may accept at least one sealed, moisture-tight, dry container of a dose of tiotropium, e.g. Spiriva®, and deliver said dose with a consistent FPD, over the expected shelf life of the product.

In a further aspect of the invention tiotropium may be mixed or formulated with at least one additional pharmacologically active ingredient(s) with an object of combining tiotropium with other medicament(s) to be used in the treatment of respiratory disorders. The present invention encompasses such use of tiotropium in a combination of medicaments directly loaded into a

sealed, moisture-tight, dry container for insertion into a DPI, the combination adapted for inhalation by the user.

The present medical product is set forth by the independent claims 1 and 2 and the dependent claims 3 to ¹²13, and a pharmaceutical combination is set forth by the independent claims ¹³14 and ¹⁴15 and the dependent claims ¹⁵16 to ²⁵24.

BRIEF DESCRIPTION OF THE DRAWINGS

The invention, together with further objects and advantages thereof, may best be understood by referring to the following detailed description taken together with the accompanying drawings, in which:

FIG. 1 illustrates in a graph the results of tests S1 to S5 and HBS1 to HBS3;

FIG. 2 illustrates in top and side views a first embodiment of a dose deposited onto a dose bed and a high barrier seal; and

FIG. 3 illustrates in top and side views a second embodiment of a dose onto a dose bed and a high barrier seal.

DETAILED DESCRIPTION

Tiotropium is a new important anticholinergic substance for treatment of asthma and COPD but tiotropium is known in the industry to have problems maintaining in-use stability due to sensitivity to moisture. This fact is also documented in the report 'COLLEGE TER BEOORDELING VAN GENEESMIDDELEN MEDICINES EVALUATION BOARD; PUBLIC ASSESSMENT REPORT; Spiriva 18 µg, inhalation powder in hard capsules; RVG 26191' (2002-05-21) on page 6/28 under 'Product development and finished product' a very short in-use stability of the Spiriva® product (9 days)

is reported and a brittleness of the capsule in the blister pack and a very low FPD: 'about 3 ug'.

Details about an inhalation kit comprising inhalable powder of tiotropium and use of an inhaler for the administration of tiotropium may also be studied in the international publication WO 03/084502 A1. Details about tiotropium compounds, medicaments based on such compounds, the use of compounds and processes for preparing compounds may be studied in the European Patent Application 0 418 716 B1.

In the light of the above information given in the quoted report a test program was set up for the physical stability of the Spiriva® product with respect to the compatibility of the formulation together with the components of the device according to Food and Drug Administration (FDA) 'Guidance for Industry; Metered Dose Inhaler (MDI) and Dry Powder Inhaler (DPI) Drug Products; Chemistry, Manufacturing, and Controls Documentation' page 37/62 'Drug product stability' lines 1209 - 1355. In 'Guidance for Industry; Stability Testing of Drug Substances and Drug Products; DRAFT GUIDANCE; B. Container/Closure' pages 35 and 36/110 lines 1127 - 1187, FDA states: 'Stability data should be developed for the drug product in each type of immediate container and closure proposed for marketing, promotion, or bulk storage. The possibility of interaction between the drug and the container and closure and the potential introduction of extractables into the drug product formulations during storage should be assessed during container/closure qualification studies using sensitive and quantitative procedures.' and further 'Loss of the active drug substance or critical excipients of the drug product by interaction with the container/closure components or components of the drug delivery device is generally evaluated as part of the stability protocol. This is usually accomplished by assaying those critical drug product components, as well as monitoring various critical parameters (e.g., pH, preservative, effectiveness). Excessive loss of a

component or change in a parameter will result in the failure of the drug product to meet applicable specifications.'

According to FDA publication 'Guidance for Industry; Stability Testing of Drug Substances and Drug Products' a 3 week test program in accelerated conditions (40 ± 2 ° / 75 ± 5 RH) for the container closure of the Spiriva® product in this case the capsule and the blister pack and the impact of the capsule and the blister package on the FPD was set up and tested.

10 Execution of tests

Spiriva® powder formulation in bulk and Spiriva® capsules from our local pharmacy where introduced to the laboratory together with the HandiHaler®. The laboratory was set up to perform in-vitro tests according to European Pharmacopoeia (EP) and US Pharmacopoeia (USP) using two Andersen cascade impactors. All analytical work where then performed according to standardized methods for Physical Tests and Determinations for Aerosols, metered-dose inhalers and dry powder inhalers described in pharmacopoeias (e.g. USP 2002 <601>) using a state of the art High Performance Liquid Chromatograph (HPLC) system.

20

Spiriva® tests

Test S1

Aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using Spiriva® formulation from bulk powder loaded into originator capsules during relative humidity below 10 %. The test was performed with 4 kPa pressure drop over the HandiHaler® at room temperature and laboratory ambient conditions.

Test S2

Aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using commercial Spiriva® capsules purchased from our local pharmacy. Test performed with 4 kPa pressure drop over

the HandiHaler® at room temperature and laboratory ambient conditions.

Test S3

5 An in-use stability test of the aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using commercial Spiriva® capsules purchased from our local pharmacy. From the blister holding 5 capsules one capsule was withdrawn and the remaining 4 capsules were put 4 days into 40 °C and 75 % Rh. The
10 blister containing the 4 capsules was then put in an excicator for 2 h before tests were performed. The test was performed with 4 kPa pressure drop over the HandiHaler® at room temperature and laboratory ambient conditions.

Test S4

15 An in-use stability test of the aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using commercial Spiriva® capsules purchased from our local pharmacy. From the blister holding 5 capsules one capsule was withdrawn and the
20 remaining 4 capsules were put 13 days into 40 °C and 75 % Rh. The blister containing the 4 capsules was then put in an excicator for 2 h before tests were performed. The test was performed with 4 kPa pressure drop over the HandiHaler® at room temperature and laboratory ambient conditions.

Test S5

25 An in-use stability test of the aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using commercial Spiriva® capsules purchased from our local pharmacy. From the
30 blister holding 5 capsules one capsule was withdrawn and the remaining 4 capsules were put 21 days into 40 °C and 75 % Rh. The blister containing the 4 capsules was then put in an excicator for 2 h

before tests were performed. The test was performed with 4 kPa pressure drop over the HandiHaler® at room temperature and laboratory ambient conditions.

5 High barrier seal tests

Test HBS1

10 An in-use stability test of the aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using Spiriva® formulation from bulk powder loaded during relative humidity below 10 % into containers made to act as a high barrier seal, in this case aluminum foils from Alcan Singen Germany and then sealed to absolute tightness. The aluminum containers were put in an exicator for 2 h before the Spiriva® powder formulation was loaded from the aluminum containers into the originator capsules at a relative
15 humidity below 10 %. The test was performed with 4 kPa pressure drop over the HandiHaler® at room temperature and laboratory ambient conditions.

Test HBS2

20 An in-use stability test of the aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using Spiriva® formulation from bulk powder loaded during relative humidity below 10 % into containers made to act as a high barrier seal, in this case aluminum foils from Alcan Singen Germany and then sealed to absolute tightness. The sealed aluminum containers were put into climate chambers for 7 days at 40 °C and 75 % Rh. The aluminum
25 containers were put in an exicator for 2 h before the Spiriva® powder formulation was loaded from the aluminum containers into the originator capsules at a relative humidity below 10 %. The test was performed with 4 kPa pressure drop over the HandiHaler® at room
30 temperature and laboratory ambient conditions.

Test HBS3

An in-use stability test of the aerodynamic fine particle fraction of metered and delivered dose out of Handihaler® using Spiriva® formulation from bulk powder loaded during relative humidity below 10 % into containers made to act as a high barrier seal, in this case aluminum foils from Alcan Singen Germany and then sealed to absolute tightness. The sealed aluminum containers were put into climate chambers for 14 days at 40 °C and 75 % Rh. The aluminum containers were then put in an exicator for 2 h before the Spiriva® powder formulation was loaded from the aluminum containers into the originator capsules at a relative humidity below 10 %. The test was performed with 4 kPa pressure drop over the HandiHaler® at room temperature and laboratory ambient conditions.

15 C-haler DPI tests

A test was also made outside the stability test program to evaluate our proprietary inhaler, the so-called C-haler, in comparison with the HandiHaler® using a tiotropium formulation. The C-haler cartridge used high barrier seals made out of aluminum foils from Alcan Singen Germany and the containers where filled volumetrically with 5 mg of the Spiriva® powder formulation in bulk. The test was performed using a 4 kPa pressure drop over the C-haler at room temperature and laboratory ambient conditions. The results from the Andersen impactor tests were calculated on fine particle fraction based on delivered dose as well as on metered dose and converted to FPD. The results are given in Table 1 below.

The results of tests S1-5 and HBS1-3 are plotted in Figure 1. The Y-axis is designated '% of commercial Spiriva® FPD'. This relates to the FPD out from the Handihaler®, where 100 % is the FPD from a fresh sample from the pharmacy.

Table 1. Inhaled fine particle dose (FPD) <5 µm in %

Calculation based on	Spiriva® in HandiHaler®, commercial sample, FPD	Spiriva® in C-haler, FPD
Metered dose	18 %	47 %
Delivered dose	36 %	56 %

5 Conclusion of the tests performed on Spiriva®

Surprisingly we have found and concluded in our tests that tiotropium is extremely sensitive to moisture and that a conventional packaging into gelatin capsules used for a majority of respiratory products will seriously affect the FPD. The results show that there is a need for a dry, moisture-tight high barrier seal enclosing the tiotropium formulation to preserve the original fine particle fraction. Not so surprisingly in the light of these findings, we have also found that the tiotropium formulation must be properly protected also during the in-use time if further reduction of the FPD shall be avoided. Eliminating the gelatin capsule has an unexpected, big, positive effect on the performance of the Spiriva® formulation.

The tests carried out show that the moisture content of the gelatin capsule reduces the FPD out of the HandiHaler® with approximately 50 % from the time of loading the dose into a capsule until the point in time when the product reaches the market. Loading Spiriva® doses into dry containers made of materials presenting high barrier seal properties and then storing the loaded containers in 40 °C and 75 % Rh, before transferring the Spiriva® doses to originator capsules and performing the same tests using HandiHaler® as before, no change can be detected in the fine particle dose (FPD), even after long periods of time. The FPD of Spiriva® in gelatin capsules, however, is further diminishing during the in-use time of the product and the FPD has been shown to drop up to another 20 % after 5 days of storage in 40 °C and 75 % Rh in an in-use stability test, due to the

breaking of the moisture barrier in the opened blister secondary package. Table 1 shows that our proprietary C-haler using high barrier containers shows a 2.6 times higher performance than HandiHaler® with respect to FPD based on metered dose.

5

State of the art

Metered doses of the Spiriva® powder formulation are today at the originator manufacturing site loaded into gelatin capsules. A gelatin capsule contains typically 13-14 % water by weight in the dose forming stage and after the capsules have been loaded they are dried in a special process in order to minimize water content. A number of dried capsules are then put in a common blister package. Details about suitable state-of-the-art capsule materials and manufacturing processes may be studied in the German Patent Application DE 101 26 924 A1. The remaining small quantity of water in the capsule material after drying is thus enclosed in the blister package and some water will be released into the enclosed air, raising the relative humidity in the air. The equilibrium between the captured air inside the package and the gelatin capsule will generate a relative humidity inside the blister package that will negatively affect the FPD of tiotropium powder out of the dry powder inhaler.

20

25

30

It is interesting to note that the big majority of dry powder formulations of many kinds of medicaments are not seriously affected by enclosed moisture in the capsule material or by normal storage variations in the relative humidity of the surrounding air. Surprisingly, our investigation has shown tiotropium to be very much different. Tiotropium powder is very much affected by very small amounts of water such that it tends to stick to wall surfaces and to agglomerate. By some mechanisms the FPD becomes less over time. Since the capsules are only used as convenient, mechanical carriers of Spiriva® doses, a solution to the moisture problem would be not to use capsules at all, but rather to directly load doses into containers made of

dry packaging material with high barrier seal properties during dry ambient conditions, preferably below 10% Rh.

The present invention discloses a dry, moisture-tight, directly loaded and sealed container enclosing a metered dose of tiotropium powder or a pharmaceutically acceptable salt, enantiomer, racemate, hydrate, or solvate, including mixtures thereof, and particularly tiotropium bromide, optionally further including excipients. The term "tiotropium" is in this document a generic term for all active forms thereof, including pharmaceutically acceptable salts, enantiomers, racemates, hydrates, solvates or mixtures thereof and may further include excipients for whatever purpose. The container uses dry, high barrier seals impervious to moisture and other foreign matters and is adapted for insertion into a dry powder inhaler device or the container may be adapted to be a part of an inhaler device.

"Dry" means that the walls of the container are constructed from selected materials such that the walls, especially the inside wall of the container, cannot release water that may affect the tiotropium powder in the dose such that the FPD is reduced. As a logical consequence container construction and materials should not be selected among those suggested in the German publication DE 101 26 924 A 1.

"High barrier seal" means a dry packaging construction or material or combinations of materials. A high barrier seal is characterized in that it represents a high barrier against moisture and that the seal itself is 'dry', i.e. it cannot give off measurable amounts of water to the load of powder. A high barrier seal may for instance be made up of one or more layers of materials, i.e. technical polymers, aluminum or other metals, glass, siliconoxides etc that together constitutes the high barrier seal.

2003 -12- 0 3

Huvudfaxen Kassa

13

A "high barrier container" is a mechanical construction made to harbour and enclose a dose of e.g. tiotropium. The high barrier container is built using high barrier seals constituting the walls of the container.

5 "Directly loaded" means that the metered dose of tiotropium is loaded directly into the high barrier container, i.e. without first loading the dose into e.g. a gelatin capsule, and then enclosing one or more of the primary containers (capsules) in a secondary package made of a high barrier seal material.

10

The high barrier containers to be loaded with tiotropium should preferably be made out of aluminum foils approved to be in direct contact with pharmaceutical products. Aluminum foils that work properly in these aspects generally consist of technical polymers laminated with aluminum foil to give the foil the correct mechanical properties to avoid cracking of the aluminum during forming. Sealing of the formed containers is normally done by using a thinner cover foil of pure aluminum or laminated aluminum and polymer. The container and cover foils are then sealed together using at least one of several possible methods, for instance:

20

- using a heat sealing lacquer, through pressure and heat;
- using heat and pressure to fuse the materials together;
- ultrasonic welding of the materials in contact.

25

Tiotropium in pure form is a very potent drug and it is therefore normally diluted before dose forming by mixing with physiologically acceptable excipients, e.g. lactose, in selected ratio(s) in order to fit a preferred method of dose forming or loading. Details about inhalation powders containing tiotropium in mixtures with excipients, methods of powder manufacture, use of powder and capsules for powder may be studied in the international publication WO 02/30389 A1.

30

In a further aspect of the invention tiotropium may be mixed or formulated with one or more other pharmacologically active ingredient(s) with an object of combining tiotropium with other medicament(s) to be used in a treatment of respiratory disorders. The present invention encompasses such use of
5 tiotropium when a combination of tiotropium and other medicaments are deposited and sealed into a dry, moisture-tight high barrier container intended for insertion into a DPI for inhalation by the user. Examples of interesting combinations of substances together with tiotropium could be inhalable steroids, nicotinamide derivatives, beta-agonists, beta-mimetics,
10 anti-histamines, adenosine A2A receptors, PDE4 inhibitors, dopamine D2 receptor agonists.

The sealed, dry, high barrier container of the invention that is directly loaded with a formulation of tiotropium may be in the form of a blister and it may
15 e.g. comprise a flat dose bed or a formed cavity in aluminum foil or a molded cavity in a polymer material, using a high barrier seal foil against ingress of moisture, e.g. of aluminum or a combination of aluminum and polymer materials. The sealed, dry, high barrier container may form a part of an inhaler device or it may be a separate item intended for insertion into an
20 inhaler device for administration of doses.

An inhaler providing a prolonged delivery of a dose during the course of a single inhalation constitutes a preferred embodiment of an inhaler for the delivery of the tiotropium powder formulation, e.g. Spiriva®. An Air-razor
25 method as described in our publication US 2003/0192539 A1 is preferably applied in the inhaler to efficiently and gradually aerosolize the dose when delivered to the user. Surprisingly enough, applying an inhaler for a prolonged delivery and using the Air-razor method on a dose comprising tiotropium in Spiriva® formulation results in a FPD at least twice as big as
30 that from the state-of-the-art HandiHaler®.

Ink. t. Patent- och req.verket

2003-12-03

Huvudfoxen Kassar

15
CLAIMS

1. A medical product comprising a dry powder dose of tiotropium directly loaded into a container, **characterized in** that

5 a dry, high barrier seal constitutes the container;
the high barrier seal of the container prevents ingress of moisture whereby the original fine particle fraction of the powder dose is preserved, and

10 the dry powder dose in the container is adapted for administration by a dry powder inhaler.

2. A medical product comprising tiotropium and separate or together with at least one additional active pharmaceutical ingredient and optionally including excipients in a dry powder medical combination dose
15 directly loaded into a container, **characterized in** that

a dry, high barrier seal constitutes the container;
the high barrier seal of the container prevents ingress of moisture whereby the original fine particle fraction of the combination dose is preserved;

20 the combination dose is adapted for administration by a dry powder inhaler, and

the at least one additional active pharmaceutical ingredient is selected from the following groups of substances: inhalable steroids, nicotinamide derivatives, beta-agonists, beta-mimetics, anti-histamines, adenosine A2A receptors, PDE4 inhibitors, dopamine D2 receptor agonists.
25

3. The medical product according to claims 1 or 2, **characterized in** that

administration of the dry powder dose is performed by inhalation
30 from a dry powder inhaler providing a prolonged dose delivery.

4. The medical product according to claims 1 or 2, **characterized in**
that

the tiotropium substance consists of one or more physiologically
acceptable salts of tiotropium.

5

5. The medical product according to claims 1 or 2, **characterized in**
that

an included excipient is lactose.

10

7. The medical product according to claims 1 or 2, **characterized in**
that

formed or flat aluminum foils, optionally laminated with
polymers, constitute the dry, high barrier seal.

15

8. The medical product according to claims 1 or 2, **characterized in**
that

a cavity molded from a polymer material selected to give the
container high barrier seal properties constitutes the container.

20

9. The medical product according to claims 1 or 2, **characterized in**
that

a cavity molded from a polymer material together with a high
barrier seal constitute the container providing it with high barrier seal
properties.

25

10. The medical product according to claims 1 or 2, **characterized in**
that

the container is a part of a dry powder inhaler.

30

11. The medical product according to claims 1 or 2, **characterized in**
that

the container is a separate part adapted for insertion into a dry powder inhaler

12. The medical product according to claims 1 or 2, **characterized in**
5 that

the container is a separate part comprising a primary part adapted for insertion into a dry powder inhaler and a secondary part enclosing the primary part in a moisture-tight package.

10 13. The medical product according to claims 1 or 2, **characterized in**
that

the medical product is intended for use in the treatment of respiratory disorders.

15 14. A pharmaceutical composition comprising tiotropium or a physiologically acceptable salt thereof and a physiologically acceptable excipient, **characterized in** that

the composition is directly loaded and sealed into a dry, moisture-tight package or dry, high barrier container in order to preserve the
20 original fine particle fraction of the composition.

25 15. A pharmaceutical composition comprising tiotropium or a physiologically acceptable salt thereof and separate or together with at least one active pharmaceutical ingredient, optionally including physiologically acceptable excipients in a combination dry powder dose, **characterized in**
that

the combination dry powder dose is directly loaded and sealed into a dry, moisture-tight package or dry, high barrier container in order to preserve the original fine particle fraction (FPF) of the composition;

30 the at least one active pharmaceutical ingredient is selected from the following groups of substances: inhalable steroids, nicotinamide

derivatives, beta antagonists, beta-mimetics, anti-histamines, adenosine A2A receptors, PDE4 inhibitors, dopamine D2 receptor agonists.

16. The pharmaceutical composition according to claim 14 or 15,
5 **characterized in that**

administration of the combination dry powder dose is performed
by inhalation from a dry powder inhaler providing a prolonged dose delivery.

17. The pharmaceutical composition according to claim 14 or 15,
10 **characterized in that**

the tiotropium comprises one or more physiologically acceptable
salts of tiotropium.

18. The pharmaceutical composition according to claim 14 or 15,
15 **characterized in that**

an included excipient is lactose.

19. The pharmaceutical composition according to claim 14 or 15,
characterized in that

20 formed or flat aluminum foils, optionally laminated with
polymers, constitute the dry, moisture-tight package or dry, high barrier
container.

20. The pharmaceutical composition according to claim 14 or 15,
25 **characterized in that**

a cavity molded from a polymer material with high barrier seal
properties forms the dry, moisture-tight package or dry, high barrier
container.

21. A pharmaceutical composition according to claim 14 or 15,
characterized in that

a cavity molded from a polymer material together with a high
barrier seal constitute the dry, moisture-tight package or dry, high barrier
5 container, thereby giving the package or container high barrier seal
properties.

22. The pharmaceutical composition according to claim 14 or 15,
characterized in that

10 the dry, moisture-tight package or dry, high barrier container
constitutes a part of a dry powder inhaler.

23. The pharmaceutical composition according to claim 14 or 15,
characterized in that

15 the dry, moisture-tight package or dry, high barrier container is a
separate part adapted for insertion into a dry powder inhaler.

24. The pharmaceutical composition according to claim 14 or 15,
characterized in that

20 the dry, moisture-tight package or dry, high barrier container is a
separate part comprising a primary package adapted for insertion into a dry
powder inhaler, and a secondary moisture-tight package or container
enclosing the primary package.

25. The pharmaceutical composition according to claim 14 or 15,
characterized in that

the pharmaceutical composition is for a use in the treatment of a
respiratory disorder.

ABSTRACT

5 The invention discloses a medical product for use in a treatment of
respiratory disorders, and comprises a metered dose of a tiotropium dry
powder formulation, directly loaded and sealed into a container made to act
as a dry high barrier seal to prevent the capture and ingress of moisture into
the tiotropium powder. The dose of tiotropium is further adapted for
10 inhalation and the container is so tight that the efficacy of the dose when
delivered is unaffected by moisture. In a further aspect of the invention a
type of inhaler is illustrated, which may accept at least one sealed, moisture-
tight container of a dose of tiotropium, to deliver the dose with a consistent
fine particle dose, over the expected shelf life of the product.

Ink. t. Patent- och reg.verket

2004-12-03

Huvudfaxen Kassa

1/2


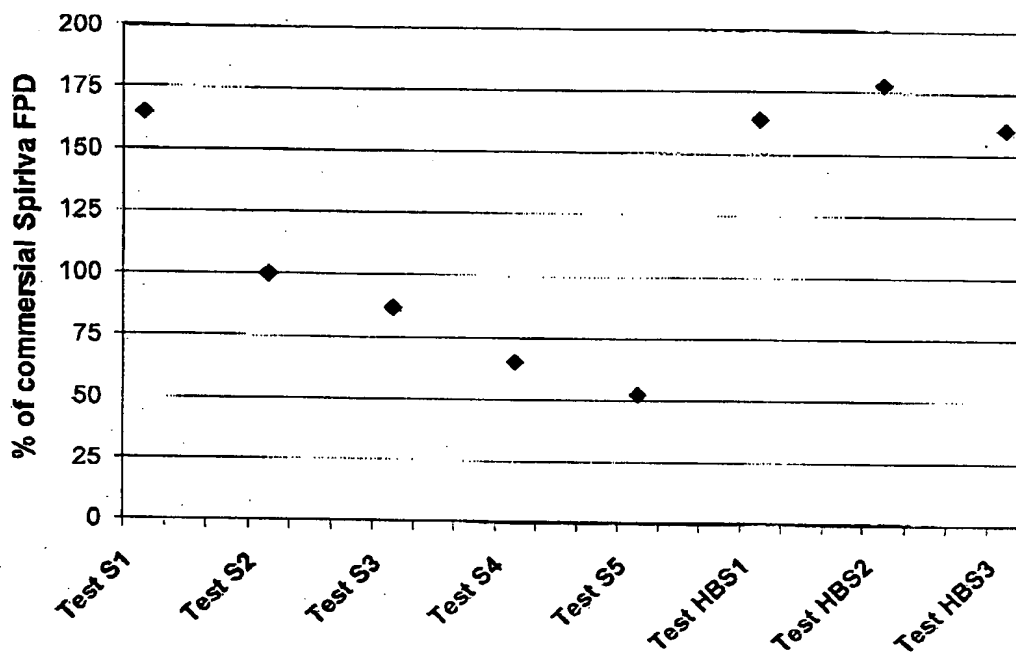


Fig. 1

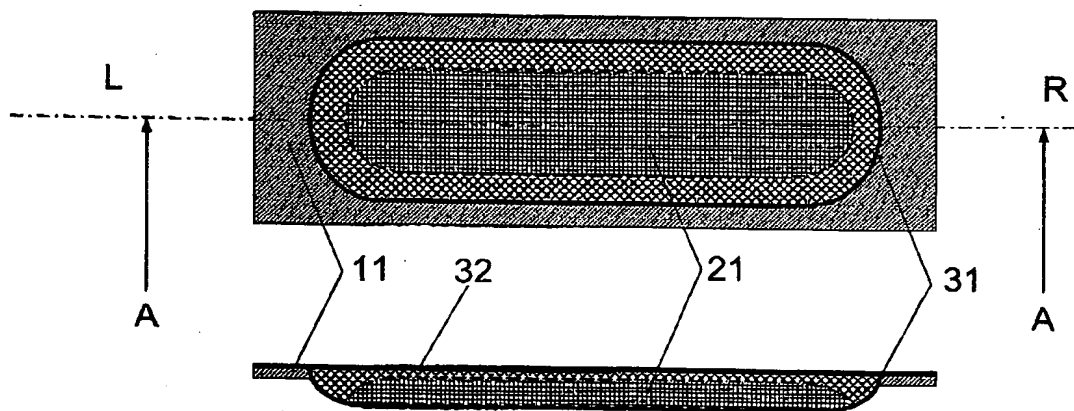

A-A

Fig. 2

2/2

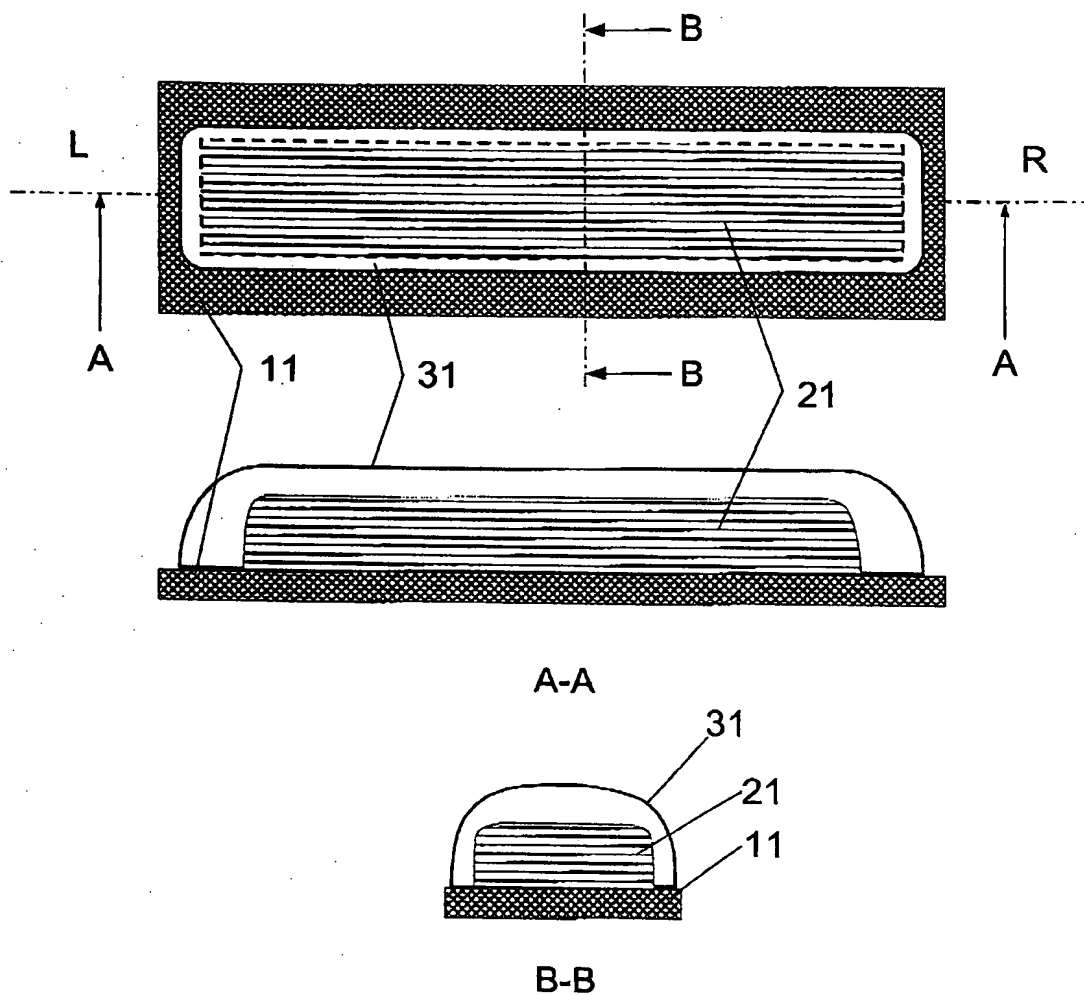


Fig. 3

BEST AVAILABLE COPY